

# BT GATT Protocol Specification - BQ Device

## Overview

This document defines the Bluetooth Low Energy (BLE) GATT services and characteristics for the BQ pill dispenser device. Uses the [BLEOTA](#) library for OTA firmware updates.

Custom services use 128-bit UUIDs with the base format: `XXXXXXXX-0000-1000-8000-00805F9B34FB`.

## Service Summary

| Service Name             | Service UUID                         | Description                              |
|--------------------------|--------------------------------------|------------------------------------------|
| Device Information (DIS) | 0x180A                               | Device identification (standard BLE DIS) |
| BLE OTA                  | 00008018-0000-1000-8000-00805f9b34fb | Firmware update via BLEOTA protocol      |
| Time Sync                | 0x1802                               | Device clock synchronization             |
| Terminal Command         | 0x1803                               | Debug/diagnostic terminal                |
| Device Sensors           | 0x1804                               | Sensor status monitoring                 |
| Device Control           | 0x1805                               | Movement and control commands            |
| Extract Pill             | 0x1806                               | Pill extraction operations               |
| Sub Drive                | 0x1807                               | Sub drive information                    |
| Notifications            | 0x1808                               | Event notifications                      |
| Offline Sync             | 0x1809                               | Offline data synchronization             |
| Battery                  | 0x180F                               | Battery level (standard BLE)             |
| Sound File               | 0x180B                               | Audio file management                    |

### 1. Device Information Service (DIS)

**Service UUID:** 0x180A (Standard BLE DIS)

Provides device identification and status information. Uses the standard Bluetooth Device Information Service combined with BLEOTA library methods.

## Characteristics

| Characteristic    | UUID   | Properties | Description              | API Method        |
|-------------------|--------|------------|--------------------------|-------------------|
| Model Number      | 0x2A24 | Read       | Device model string      | setModel()        |
| Serial Number     | 0x2A25 | Read       | Device serial number     | setSerialNumber() |
| Firmware Revision | 0x2A26 | Read       | Firmware version string  | setFWVersion()    |
| Hardware Revision | 0x2A27 | Read       | Hardware revision string | setHWVersion()    |
| Manufacturer Name | 0x2A29 | Read       | Manufacturer identifier  | setManufactuer()  |

## Battery Service

Battery level is provided via a separate standard Battery Service ( 0x180F ):

| Characteristic | UUID   | Properties   | Description                        |
|----------------|--------|--------------|------------------------------------|
| Battery Level  | 0x2A19 | Read, Notify | Current battery percentage (0-100) |

## Battery Level Format

Byte 0: Battery percentage (0x00 - 0x64)

## Firmware Version Format

Bytes 0-N: Major.Minor.Patch as ASCII string (null-terminated)  
Example: "1.2.3"

---

## 2. BLE OTA Service (BLEOTA)

**Service UUID:** 00008018-0000-1000-8000-00805f9b34fb

**Reference:** [github.com/gb88/BLEOTA](https://github.com/gb88/BLEOTA)

Handles over-the-air (OTA) firmware and SPIFFS filesystem updates using the BLEOTA protocol. Supports zlib compression (~30% faster transfers) and RSA-2048 signature

verification.

## Characteristics

| Characteristic   | UUID                                 | Properties             | Description                      |
|------------------|--------------------------------------|------------------------|----------------------------------|
| Firmware Receive | 00008020-0000-1000-8000-00805f9b34fb | Write, WriteNR, Notify | Firmware data reception with ACK |
| Command          | 00008022-0000-1000-8000-00805f9b34fb | Write, WriteNR, Notify | OTA commands and responses       |

## Command Packet Format

Bytes 0-1: Command ID (uint16, little-endian)

Bytes 2-17: Payload (16 bytes, zero-padded if unused)

Bytes 18-19: CRC16 checksum

Total: 20 bytes

## Command Codes

| Command          | ID     | Payload                             | Description                 |
|------------------|--------|-------------------------------------|-----------------------------|
| Start Flash OTA  | 0x0001 | Bytes 2-5: Firmware size (uint32)   | Initialize firmware update  |
| Start SPIFFS OTA | 0x0004 | Bytes 2-5: Filesystem size (uint32) | Initialize SPIFFS update    |
| Stop OTA         | 0x0002 | None                                | Abort current update        |
| Response         | 0x0003 | Status code                         | Device response to commands |

## Command Response Status Codes

| Status          | Value  | Description                       |
|-----------------|--------|-----------------------------------|
| Accept          | 0x0000 | Command accepted, OTA started     |
| Reject          | 0x0001 | Command rejected                  |
| Signature Error | 0x0003 | RSA signature verification failed |

## Firmware Data Packet Format (App → Device)

Bytes 0-1: Sector Index (uint16, little-endian)  
Byte 2: Packet Sequence (0x00-0xFE for data, 0xFF for final packet)  
Bytes 3-N: Payload data (MTU - 4 bytes, typically up to 243 bytes)

Note: When Packet\_Seq = 0xFF, the final 2 bytes contain CRC16 of the sector

Firmware ACK Format (Device → App via Notify)


Bytes 0-1: Sector Index (uint16, little-endian)  
Bytes 2-3: ACK Status (uint16, little-endian)  
Bytes 4-5: CRC16 checksum


Total: 6 bytes

ACK Status Codes

| Status            | Value  | Description                  |
|-------------------|--------|------------------------------|
| Success           | 0x0000 | Sector received successfully |
| CRC Error         | 0x0001 | CRC verification failed      |
| Sector Mismatch   | 0x0002 | Unexpected sector index      |
| Length Error      | 0x0003 | Invalid payload length       |
| OTA Start Failure | 0x0005 | Failed to initialize OTA     |

OTA Update Flow

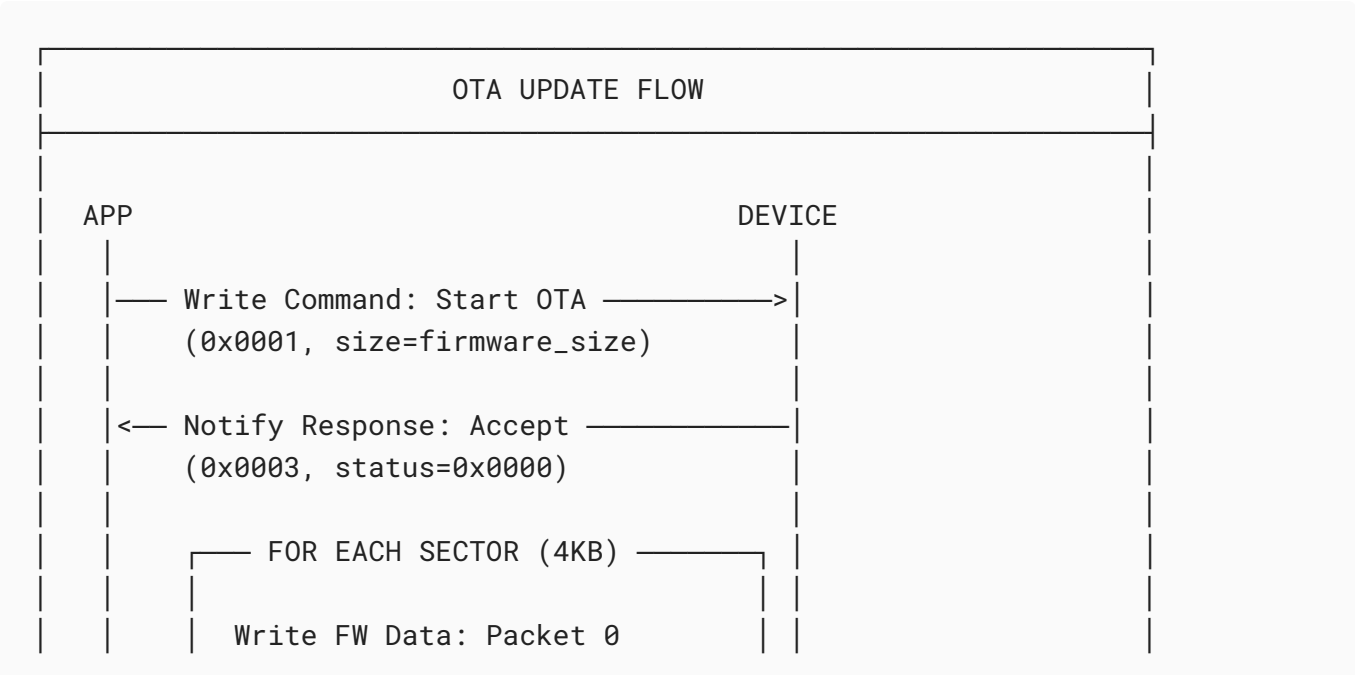

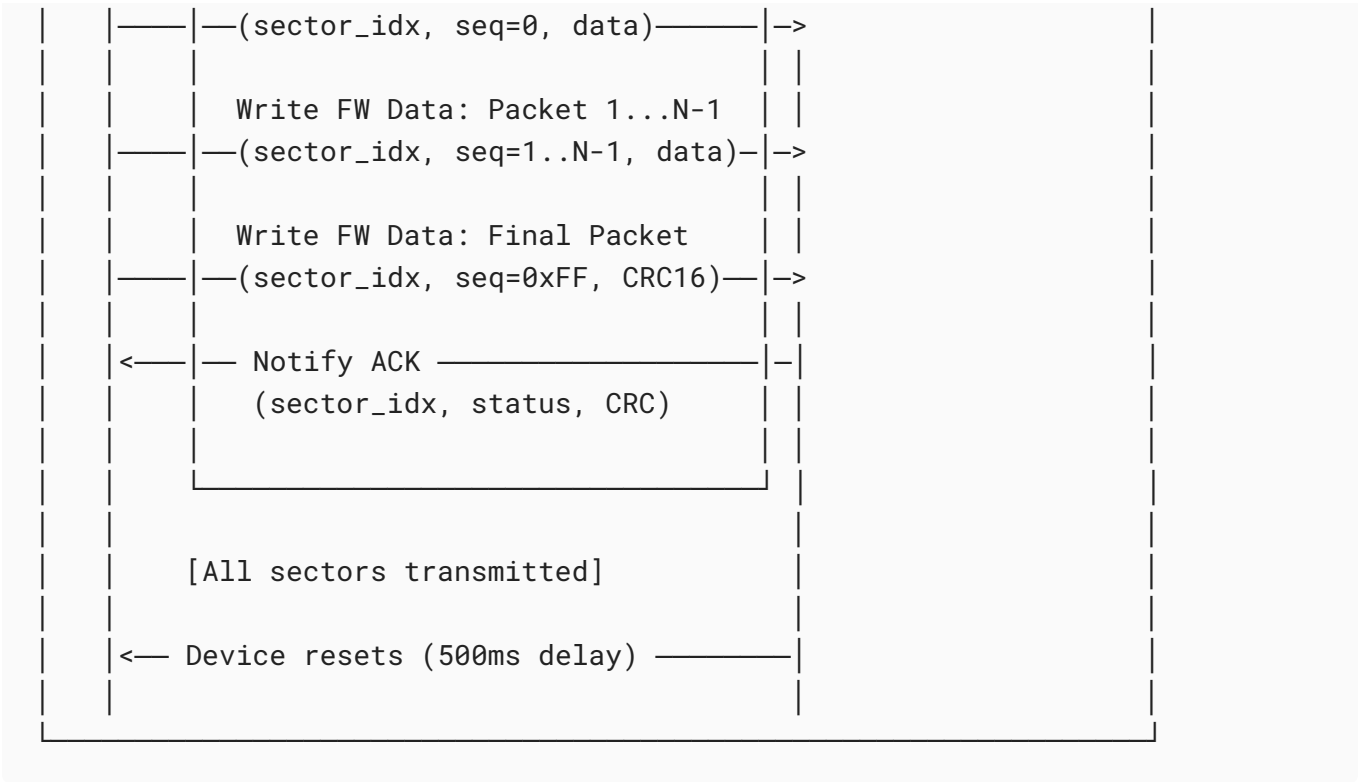

## Sector Transmission Details

| Parameter          | Value                                            |
|--------------------|--------------------------------------------------|
| Sector Size        | 4096 bytes (4KB)                                 |
| Packet Size        | MTU - 4 bytes (typically 243 bytes with MTU=247) |
| Packets per Sector | ~17 packets for 4KB sector                       |
| Sequence Numbers   | 0x00 to 0xFE for data, 0xFF signals final packet |

## Security Features

| Feature     | Description                                                      |
|-------------|------------------------------------------------------------------|
| CRC16       | Per-sector integrity check                                       |
| RSA-2048    | Firmware signature verification using SHA256 hash                |
| Compression | Optional zlib compression for faster transfer (~30% improvement) |

## Post-Update Behavior

After successful OTA completion:

1. Device validates firmware signature (if enabled)

2. Device waits 500ms (configurable)
3. Device resets and boots new firmware

## Error Recovery

| Scenario        | Recovery                           |
|-----------------|------------------------------------|
| CRC Error       | Retransmit current sector          |
| Timeout         | Restart from last ACK'd sector     |
| Signature Fail  | OTA rejected, device unchanged     |
| Connection Lost | Resume from last successful sector |

## Example: Calculate Transfer Time

```
Firmware size: 1MB (1,048,576 bytes)
Sector size: 4KB
Total sectors: 256
Packets per sector: ~17
Connection interval: 30ms
Packets per interval: 6 (conservative)

Time per sector: 17 / 6 * 30ms ≈ 85ms
Total time: 256 * 85ms ≈ 22 seconds

With zlib compression (~30% reduction):
Compressed size: ~734KB
Total time: ≈ 15 seconds
```

---

## 3. Time Sync Service

**Service UUID:** 0x1802

Synchronizes device clock with mobile app/server time.

### Characteristics

| Characteristic | UUID   | Properties  | Description         |
|----------------|--------|-------------|---------------------|
| Current Time   | 0x2C01 | Read, Write | Device current time |
| Time Zone      | 0x2C02 | Read, Write | Timezone offset     |

| Characteristic | UUID   | Properties  | Description            |
|----------------|--------|-------------|------------------------|
| DST Offset     | 0x2C03 | Read, Write | Daylight saving offset |
| Last Sync      | 0x2C04 | Read        | Timestamp of last sync |

## Current Time Format

Bytes 0-3: Unix timestamp (uint32, little-endian, seconds since 1970)  
 Bytes 4-5: Milliseconds (uint16, little-endian)

## Time Zone Format

Byte 0: Timezone offset in 15-minute increments (int8)  
 Example: UTC+2 = 0x08 (8 \* 15 = 120 minutes)

## DST Offset Format

Byte 0: DST offset in minutes (uint8)  
 0x00 = No DST  
 0x3C = 60 minutes DST

# 4. Terminal Command Service

**Service UUID:** 0x1803

Provides debug and diagnostic terminal access.

## Characteristics

| Characteristic  | UUID   | Properties   | Description                 |
|-----------------|--------|--------------|-----------------------------|
| TX (Device→App) | 0x2D01 | Notify       | Terminal output from device |
| RX (App→Device) | 0x2D02 | Write        | Terminal input to device    |
| Command Status  | 0x2D03 | Read, Notify | Command execution status    |

## Command Format

Bytes 0-N: ASCII command string (null-terminated, max 200 bytes)

# Response Format

```
Byte 0: Response type
  0x00 = Output line
  0x01 = Command complete (success)
  0x02 = Command complete (error)
  0x03 = Command in progress
Bytes 1-N: ASCII response data
```

## Common Terminal Commands

| Command     | Description                           |
|-------------|---------------------------------------|
| help        | List available commands               |
| status      | Device status summary                 |
| reset       | Soft reset device                     |
| factory     | Factory reset (requires confirmation) |
| log [level] | Set/get log level                     |
| diag        | Run diagnostics                       |

## 5. Device Sensors Service

**Service UUID:** 0x1804

Monitors device sensor status and readings.

### Characteristics

| Characteristic | UUID   | Properties   | Description                   |
|----------------|--------|--------------|-------------------------------|
| Sensor Status  | 0x2E01 | Read, Notify | Combined sensor status bitmap |
| Temperature    | 0x2E02 | Read, Notify | Internal temperature          |
| Humidity       | 0x2E03 | Read, Notify | Internal humidity             |
| Motor Status   | 0x2E04 | Read, Notify | Motor position and status     |
| Lid Status     | 0x2E05 | Read, Notify | Lid open/closed status        |
| Pill Detection | 0x2E06 | Read, Notify | Pill presence sensors         |



## Sensor Status Bitmap

Byte 0: Status flags

- Bit 0: Temperature sensor OK
- Bit 1: Humidity sensor OK
- Bit 2: Motor encoder OK
- Bit 3: Lid sensor OK
- Bit 4: Pill sensors OK
- Bit 5: Battery sensor OK
- Bit 6-7: Reserved

Byte 1: Error flags

- Bit 0: Temperature out of range
- Bit 1: Humidity out of range
- Bit 2: Motor fault
- Bit 3: Lid jam detected
- Bit 4: Pill sensor fault
- Bit 5-7: Reserved

## Temperature Format

Bytes 0-1: Temperature in 0.01°C (int16, little-endian)  
Example: 2350 = 23.50°C

## Motor Status Format

Byte 0: Current position (0-11 for 12 slots)

Byte 1: Target position

Byte 2: Motor state

- 0x00 = Idle
- 0x01 = Moving
- 0x02 = Homing
- 0x03 = Error

Byte 3: Error code (if state = 0x03)

## Pill Detection Format

Bytes 0-1: Pill presence bitmap (uint16)

- Bit N = 1 if pill present in slot N

Example: 0xFFFF = all 12 slots have pills

---

## 6. Device Control Service

**Service UUID:** 0x1805

Controls device movement and operations.

### Characteristics

| Characteristic  | UUID   | Properties          | Description              |
|-----------------|--------|---------------------|--------------------------|
| Control Command | 0x2F01 | Write, Indicate     | Issue control commands   |
| Control Status  | 0x2F02 | Read, Notify        | Current operation status |
| Position        | 0x2F03 | Read, Write, Notify | Current/target position  |

### Control Commands

| Command   | Value | Parameters        | Description              |
|-----------|-------|-------------------|--------------------------|
| Home      | 0x01  | None              | Return to home position  |
| Rotate    | 0x02  | Position (1 byte) | Rotate to position 0-11  |
| Extract   | 0x03  | Position (1 byte) | Extract pill at position |
| Stop      | 0x04  | None              | Emergency stop           |
| Calibrate | 0x05  | None              | Run calibration sequence |

### Command Request Format

Byte 0: Command code  
Bytes 1-N: Command parameters (varies by command)

### Command Response Format (Indicate)

Byte 0: Original command code  
Byte 1: Result code  
    0x00 = Success  
    0x01 = Invalid command  
    0x02 = Invalid parameter  
    0x03 = Busy  
    0x04 = Hardware error  
    0x05 = Timeout  
Bytes 2-N: Result data (varies by command)

# Control Status Format

```
Byte 0: Device state
    0x00 = Idle
    0x01 = Executing command
    0x02 = Error state
    0x03 = Calibrating
Byte 1: Last command code
Byte 2: Last command result
Byte 3: Current position (0-11)
```

## 7. Extract Pill Service

**Service UUID:** 0x1806

Manages pill extraction operations.

### Characteristics

| Characteristic  | UUID   | Properties      | Description                 |
|-----------------|--------|-----------------|-----------------------------|
| Extract Command | 0x3001 | Write, Indicate | Initiate extraction         |
| Extract Status  | 0x3002 | Read, Notify    | Extraction operation status |
| Extract History | 0x3003 | Read            | Recent extraction records   |
| Schedule        | 0x3004 | Read, Write     | Scheduled extractions       |

### Extract Command Format

```
Byte 0: Command type
    0x01 = Extract from position
    0x02 = Extract next scheduled
    0x03 = Cancel extraction
Byte 1: Position (0-11, for type 0x01)
Byte 2: Confirmation flag (0x5A = confirmed)
```

### Extract Status Format

```
Byte 0: Extraction state
    0x00 = Idle
```

```
0x01 = Rotating to position
0x02 = Extracting
0x03 = Verifying
0x04 = Complete (success)
0x05 = Complete (failed)
Byte 1: Current position
Byte 2: Target position
Byte 3: Error code (if failed)
    0x00 = No error
    0x01 = No pill detected
    0x02 = Jam detected
    0x03 = Motor error
    0x04 = Timeout
```

## Extract History Format

```
Each record (8 bytes):
    Bytes 0-3: Timestamp (uint32, Unix time)
    Byte 4: Position extracted
    Byte 5: Result (0x00 = success, other = error code)
    Bytes 6-7: Reserved

Max 10 records returned
```

## Schedule Format

```
Each schedule entry (6 bytes):
    Byte 0: Day mask (bit 0 = Sunday, bit 6 = Saturday)
    Byte 1: Hour (0-23)
    Byte 2: Minute (0-59)
    Byte 3: Position (0-11)
    Byte 4: Enabled flag (0x00 = disabled, 0x01 = enabled)
    Byte 5: Reserved

Max 14 schedule entries
```

---

## 8. Sub Drive Service

**Service UUID:** 0x1807

Provides sub drive (cartridge) information.

# Characteristics

| Characteristic   | UUID   | Properties   | Description              |
|------------------|--------|--------------|--------------------------|
| Sub Drive Info   | 0x3101 | Read, Notify | Cartridge identification |
| Sub Drive Status | 0x3102 | Read, Notify | Cartridge status         |
| Pill Count       | 0x3103 | Read, Notify | Pills remaining per slot |
| Medication Info  | 0x3104 | Read         | Medication metadata      |

## Sub Drive Info Format

```
Bytes 0-15: Cartridge UUID (128-bit)
Bytes 16-19: Manufacturing date (uint32, Unix time)
Bytes 20-23: Expiration date (uint32, Unix time)
Byte 24: Total slots (typically 12)
Byte 25: Cartridge type
    0x01 = Standard
    0x02 = Extended
    0x03 = Custom
```

## Sub Drive Status Format

```
Byte 0: Status
    0x00 = Not inserted
    0x01 = Inserted, valid
    0x02 = Inserted, expired
    0x03 = Inserted, error
Byte 1: Error code (if status = 0x03)
Bytes 2-3: Total extractions (uint16)
Bytes 4-5: Successful extractions (uint16)
```

## Pill Count Format

```
Bytes 0-11: Pill count per slot (1 byte each, 0-255)
    Slot 0 at byte 0, slot 11 at byte 11
Byte 12: Total pills remaining
```

## Medication Info Format

Bytes 0-31: Medication name (UTF-8, null-terminated)  
Bytes 32-47: NDC code (ASCII, null-terminated)  
Bytes 48-51: Dosage (uint32, micrograms)  
Byte 52: Dosage unit  
    0x01 = mg  
    0x02 = mcg  
    0x03 = IU

## 9. Notifications Service

**Service UUID:** 0x1808

Handles event notifications and alerts.

### Characteristics

| Characteristic      | UUID   | Properties   | Description               |
|---------------------|--------|--------------|---------------------------|
| Notification        | 0x3201 | Notify       | Push notifications to app |
| Notification Config | 0x3202 | Read, Write  | Notification settings     |
| Alert Status        | 0x3203 | Read, Notify | Active alert status       |
| Acknowledge         | 0x3204 | Write        | Acknowledge alerts        |

### Notification Format

Byte 0: Notification type  
    0x01 = Scheduled extraction reminder  
    0x02 = Extraction complete  
    0x03 = Extraction failed  
    0x04 = Low battery  
    0x05 = Cartridge empty  
    0x06 = Cartridge expiring  
    0x07 = Device error  
    0x08 = Connectivity lost  
Bytes 1-4: Timestamp (uint32, Unix time)  
Byte 5: Priority (0x01 = low, 0x02 = normal, 0x03 = high)  
Bytes 6-N: Payload (varies by type, max 100 bytes)

### Notification Config Format

Byte 0: Enabled notifications bitmap  
  Bit 0: Scheduled reminders  
  Bit 1: Extraction events  
  Bit 2: Battery alerts  
  Bit 3: Cartridge alerts  
  Bit 4: Error alerts  
  Bit 5-7: Reserved  
Byte 1: Reminder lead time (minutes before scheduled)  
Byte 2: Low battery threshold (percentage)  
Byte 3: Quiet hours start (hour, 0-23)  
Byte 4: Quiet hours end (hour, 0-23)

## Alert Status Format

Byte 0: Active alerts bitmap  
  Bit 0: Low battery  
  Bit 1: Empty slot  
  Bit 2: Expired cartridge  
  Bit 3: Motor error  
  Bit 4: Sensor error  
  Bit 5: Missed dose  
  Bit 6-7: Reserved  
Bytes 1-4: Oldest alert timestamp

## Acknowledge Format

Byte 0: Alert type to acknowledge (from alert bitmap)  
  0xFF = Acknowledge all

---

# 10. Offline Sync Service

**Service UUID:** 0x1809

Synchronizes offline data (sub drive logs, telemetry).

## Characteristics

| Characteristic | UUID   | Properties      | Description             |
|----------------|--------|-----------------|-------------------------|
| Sync Control   | 0x3301 | Write, Indicate | Control sync operations |

| Characteristic | UUID   | Properties   | Description                    |
|----------------|--------|--------------|--------------------------------|
| Sync Data      | 0x3302 | Notify       | Streamed sync data             |
| Sync Status    | 0x3303 | Read, Notify | Sync progress status           |
| Last Sync Time | 0x3304 | Read         | Last successful sync timestamp |

## Sync Control Commands

| Command       | Value | Parameters                | Description               |
|---------------|-------|---------------------------|---------------------------|
| Start Sync    | 0x01  | Start timestamp (4 bytes) | Begin sync from timestamp |
| Pause Sync    | 0x02  | None                      | Pause current sync        |
| Resume Sync   | 0x03  | None                      | Resume paused sync        |
| Abort Sync    | 0x04  | None                      | Cancel sync operation     |
| Mark Complete | 0x05  | End timestamp (4 bytes)   | Mark sync complete        |

## Sync Data Format

```
Byte 0: Record type
    0x01 = Extraction event
    0x02 = Error event
    0x03 = Telemetry sample
    0x04 = State change
Bytes 1-4: Timestamp (uint32, Unix time)
Byte 5: Sequence number (for ordering)
Bytes 6-N: Record payload (varies by type)
```

## Extraction Event Payload

```
Byte 0: Position
Byte 1: Result (0x00 = success, other = error)
Byte 2: Duration (seconds)
Byte 3: Battery level at time
```

## Telemetry Sample Payload

```
Bytes 0-1: Temperature (int16, 0.01°C)
Bytes 2-3: Humidity (uint16, 0.01%)
```



Byte 4: Battery level  
Byte 5: Motor status

## Sync Status Format

Byte 0: Sync state  
    0x00 = Idle  
    0x01 = Syncing  
    0x02 = Paused  
    0x03 = Error  
Bytes 1-2: Records pending (uint16)  
Bytes 3-4: Records synced (uint16)  
Byte 5: Progress percentage (0-100)

# 11. Sound File Service

**Service UUID:** 0x180B

Manages audio files for alerts and notifications.

## Characteristics

| Characteristic | UUID   | Properties             | Description            |
|----------------|--------|------------------------|------------------------|
| Sound Control  | 0x3401 | Write, Indicate        | Sound playback control |
| Sound Upload   | 0x3402 | Write Without Response | Upload sound data      |
| Sound List     | 0x3403 | Read                   | List available sounds  |
| Sound Config   | 0x3404 | Read, Write            | Sound configuration    |
| Volume         | 0x3405 | Read, Write            | Master volume level    |

## Sound Control Commands

| Command      | Value | Parameters               | Description           |
|--------------|-------|--------------------------|-----------------------|
| Play         | 0x01  | Sound ID (1 byte)        | Play sound by ID      |
| Stop         | 0x02  | None                     | Stop current playback |
| Upload Start | 0x03  | Sound ID, Size (4 bytes) | Begin upload          |
| Upload End   | 0x04  | Checksum (4 bytes)       | Complete upload       |

| Command | Value | Parameters        | Description       |
|---------|-------|-------------------|-------------------|
| Delete  | 0x05  | Sound ID (1 byte) | Delete sound file |

## Sound Upload Format

Bytes 0-1: Packet sequence (uint16)  
 Bytes 2-N: PCM audio data (8kHz, 8-bit mono, max 244 bytes)

## Sound List Format

Each entry (20 bytes):  
 Byte 0: Sound ID  
 Bytes 1-16: Sound name (UTF-8, null-terminated)  
 Bytes 17-18: Duration (uint16, milliseconds)  
 Byte 19: Type  
     0x01 = System  
     0x02 = Alert  
     0x03 = Custom

## Sound Config Format

Byte 0: Reminder sound ID  
 Byte 1: Success sound ID  
 Byte 2: Error sound ID  
 Byte 3: Low battery sound ID  
 Byte 4: Extraction sound ID

## Volume Format

Byte 0: Volume level (0-100)  
 Byte 1: Mute flag (0x00 = not muted, 0x01 = muted)

---

## Error Codes Reference

### Global Error Codes

| Code | Name              | Description                          |
|------|-------------------|--------------------------------------|
| 0x00 | SUCCESS           | Operation completed successfully     |
| 0x01 | ERR_INVALID_PARAM | Invalid parameter value              |
| 0x02 | ERR_INVALID_STATE | Operation not valid in current state |
| 0x03 | ERR_BUSY          | Device busy with another operation   |
| 0x04 | ERR_TIMEOUT       | Operation timed out                  |
| 0x05 | ERR_NOT_FOUND     | Requested resource not found         |
| 0x06 | ERR_NO_MEMORY     | Insufficient memory                  |
| 0x07 | ERR_HARDWARE      | Hardware fault detected              |
| 0x08 | ERR_COMM          | Communication error                  |
| 0x09 | ERR_AUTH          | Authentication required/failed       |
| 0x0A | ERR_NOT_SUPPORTED | Feature not supported                |
| 0xFF | ERR_UNKNOWN       | Unknown error                        |

---

## Connection Parameters

### Recommended BLE Parameters

| Parameter           | Value     | Description                       |
|---------------------|-----------|-----------------------------------|
| Connection Interval | 30-50 ms  | Balance between latency and power |
| Slave Latency       | 0-4       | Allow skipping connection events  |
| Supervision Timeout | 4000 ms   | Connection timeout                |
| MTU                 | 247 bytes | Maximum transmission unit         |

### Advertising Parameters

| Parameter            | Value                                  |
|----------------------|----------------------------------------|
| Advertising Interval | 100-200 ms (fast), 1000-2000 ms (slow) |
| Advertising Type     | Connectable, Scannable                 |
| TX Power             | 0 dBm                                  |

# Service Discovery

The device advertises the following service UUIDs:

- Device Information Service ( 0x180A )
- BLE OTA Service ( 00008018-0000-1000-8000-00805f9b34fb )
- Device Control Service ( 0x1805 )

Full service list available after connection via GATT discovery.

---

## Security

### Pairing Requirements

- **Bonding:** Required
- **MITM Protection:** Required
- **Secure Connections:** Supported (LE Secure Connections)
- **IO Capability:** Display + Yes/No

### Encryption

All characteristics except Device Information Service (read-only info) require an encrypted connection.

### Authentication

Sensitive operations (FW Update, Factory Reset, Extract) require:

1. Encrypted connection
  2. Bonded device
  3. Optional: PIN confirmation for critical operations
- 

## Version History

| Version | Date | Changes               |
|---------|------|-----------------------|
| 1.0     | TBD  | Initial specification |